




## STELLA, THE GASTRIC BALLOON THAT ENSURES FAST, EASY AND SAFE PLACEMENT

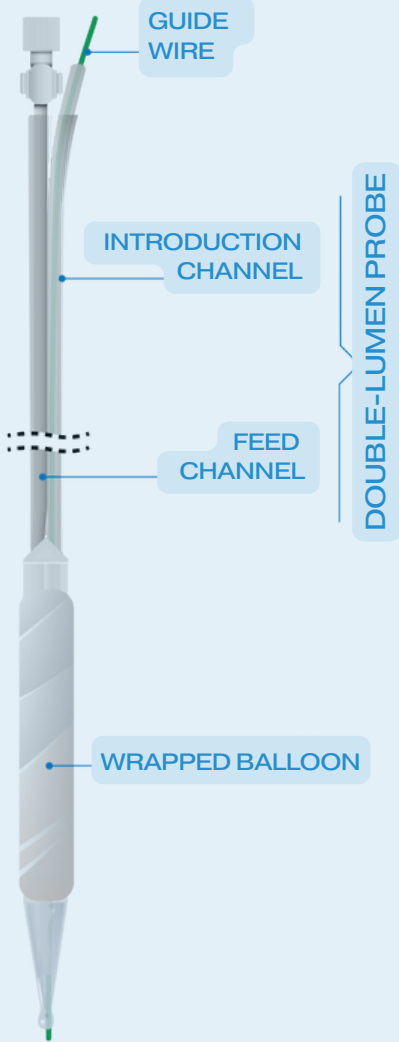

The balloon is wrapped around itself in a sleeve that opens when filled, and is fixed to the distal end of the feed channel. The sleeve, which **protects the balloon and facilitates its placement**, is attached to the introduction channel, through which the guide wire passes.

The probe that enables the balloon to be placed into the patient's gastric cavity is a double-lumen probe. One of the lumens acts as a feed channel for the balloon once the introducer has been removed, which should be placed in the balloon insertion procedure. The other lumen allows the guide wire, previously inserted by endoscopy, to slide inside it and guide the Stella device into the patient **smoothly, quickly and without misjudgement**.



Örnek Mah. Ercüment Batanay Sokak  
Dumankaya İkon Sitesi A1 Blok No:14  
Kat:1 Daire: 6 Ataşehir-İstanbul / TÜRKİYE  
+90 216 302 03 20 / [info@mendomed.com.tr](mailto:info@mendomed.com.tr)  
[www.mendomed.com.tr](http://www.mendomed.com.tr)

